


Aromatics (Benzene/Toluene/Xylene) Rtx[®]-502.2

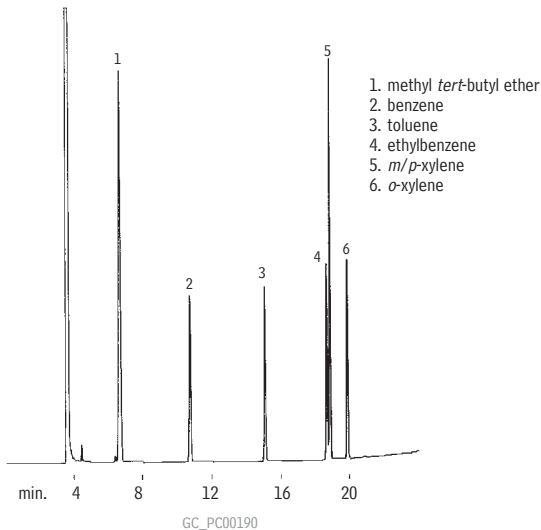

105m, 0.53mm ID, 3.0 μ m Rtx[®]-502.2 (cat.# 10910)
0.25 μ L split injection of CA PVOC Mix (cat.# 30071)
Conc.: 1000–5000 μ g/mL

Oven temp.: 40°C (hold 1 min.) to 100°C @ 5°C/min.,
to 240°C @ 8°C/min. (hold 8 min.)

Inj. / det. temp.: 200°C / 280°C
Carrier gas: helium
Linear velocity: 67cm/sec. set @ 40°C
FID sensitivity: 2.05 x 10⁹ AFS
Split flow: 70mL/min.